**NECK AND SHOULDER MASSAGER
MASAJEADOR**

**RECHARGEABLE NECK AND
SHOULDER MASSAGER
MASAJEADOR RECARGABLE**

® , WESTINGHOUSE, and INNOVATION YOU CAN BE SURE OF are trademarks of Westinghouse Electric Corporation. Used under license by Westinghouse Homeware (HK) Co. All Rights Reserved.

WHMSG2200 / WHMSGR2200 • 09/20


Instructions

For household use only

Instrucciones

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English

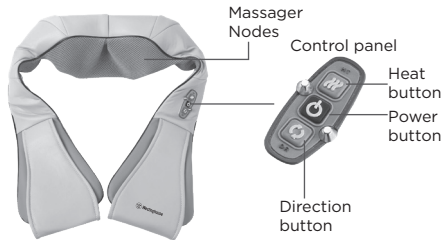INTENDED USE

Use the device only for the purpose intended. Any other use may result in damage to the device or injuries.

If you have the following disorder or condition, please consult your physician before use:

1. Heart failure, kidney failure, infection of the superficial veins (called phlebitis) or soft tissue (called cellulitis) in the legs or elsewhere, blood clots in the legs and bleeding disorders.
2. You are pregnant.
3. You have a pacemaker.
4. You have any concerns regarding your health.

IMPORTANT SAFEGUARDS

Specific safety instructions

1. Do not operate with wet hands.
2. Do not operate in water.
3. Do not drag the power jack cable.
4. Do not place heavy objects on the massage pillow.
5. Do not sleep while using the massager.
6. Do not use while intoxicated.
7. Do not directly massage hand, elbow parts, knee, breast or abdominal area.
8. Do not exceed 30 minutes of use.
9. To prevent electrical shock, DO NOT pull manual power jack cable. Disconnect power cord by the plug. Do not pull at the plug body.
10. If massager does not function, turn power off and discontinue use.
11. Do not allow children to use or play with the massager.
12. Discontinue use in the presence of discomfort.
13. Seek medical advice prior to use if you are pregnant, have a fever, malignant tumor, have heart disease, recovering from a surgical procedure, or have an acute disorder.
14. Keep away from sharp items. Repair if product's material is cut or torn.
15. Disconnect from power source before wiping with damp cloth.
16. Do not use if power cord is frayed or damaged.
17. This product does not include a charger, if you use one it must be 5V / 1A.

BEFORE FIRST USE

1. Carefully take the device and its accessories out of the packaging. Retain the packaging in case you would like to store the device away or transport it at a later date. If you choose to dispose of the packaging, then please do so according to the applicable legal regulations.
2. Check the packaging contents for completeness and damages. Should the contents in the packaging not be complete or be damaged, do not use the device. Return it to your retail outlet immediately.
3. Clean the device and all accessories before the first use. Refer to the chapter CARE AND MAINTENANCE.

USING THE NECK & SHOULDER MASSAGER.

1. press the power button to turn the massager on.
2. Press direction button to change the clockwise and anticlockwise.
3. Press heating button shortly to on or off the heating.
4. After massage press the power button to turn the massager off.

CARE AND MAINTENANCE

1. Unplug the massager and let the massager nodes cool down if just used.
2. Clean the outside of massager with a damp cloth and dry with a soft, clean cloth.
3. Avoid wiping down the mesh covering the massager nodes, this material is very thin and water will seep through it.
4. Let the massager's cover dry completely before using and make sure that the AC adaptor is not near water when using again.
5. The massager should be stored in a dry, cool location, away from sources of water or heat. It is best stored in a secure box or in an airtight bag to keep it free of dust.

DISPOSAL

Proper disposal of this product.


This label indicates that, in the EU, electronic devices cannot be disposed of with regular household waste.

Please use the return or collections systems of your municipality or contact the seller where you purchased this product.

INFORMATION AND SERVICE

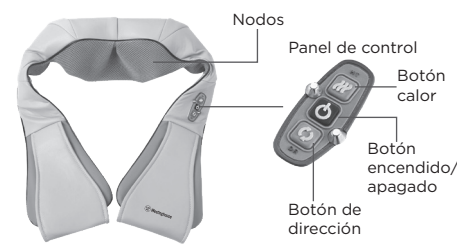If you have question or concerns about your product, please visit our website: www.westinghousehomeware.com, or contact our oversea offices/agents. In line with our policy of continuous product development, we reserve the right to change product, packaging and documentation specifications without notice.

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